



AGA - 1 shpk

Official Distributor

WELL D-ONE® SYSTEM


UROGEN WELL D-ONE®

FOR THE MOST COMMON UROGENITAL INFECTIONS

System composed of a polypropylene plate containing 32 conical wells for better viewing of the colorimetric reactions that occur as result of the growth of microorganisms in specially formulated media for selective culture of: *Mycoplasma hominis*, *Mycoplasma spp.*, *Ureaplasma urealyticum/parvum*.

The Kit allows the presumptive identification of *Neisseria spp.*, *Streptococcus B*, *Staphylococcus aureus*,

Enterococcus spp., *Escherichia coli*, *Gardnerella vaginalis*, *Trichomonas vaginalis* and *Candida spp.*

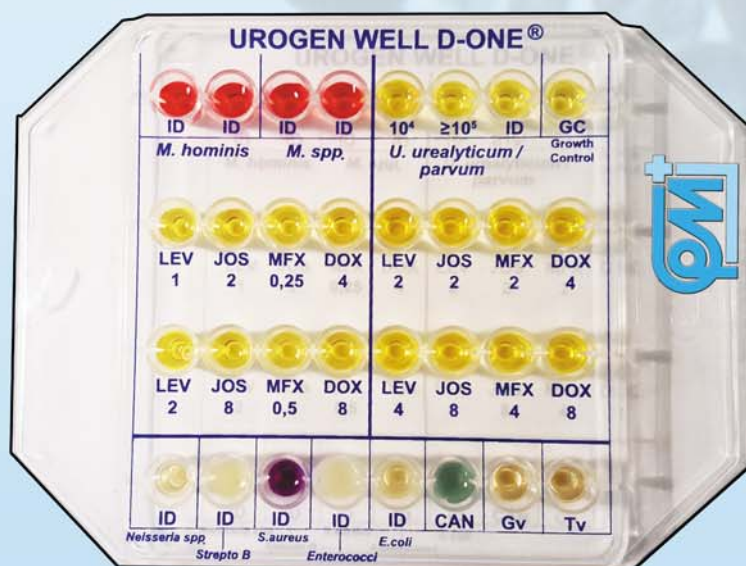

The tested sample is diluted in accordance with the procedure and inoculated directly into the wells of the plate.

Incubation is performed for 18-24-48 hours and readings are taken visually without additional equipment.



The results are evidenced by color changes that occur in wells according to reactions due to chemical or chromogenic components contained in specific formulations for each microorganism.